

# MedMinder Pill Dispenser's Features at a glance



## Reminders

Users are prompted to take or refill their medications through a series of customizable visual and auditory alerts. Set up the reminders to blink, beep or play a custom audio message. If the user doesn't respond to these reminders they will receive a phone call. If there is no response the family or caregiver are notified by phone, email or text message ensuring that everyone stays informed.

## Remote Programming And Monitoring


MedMinder is equipped with a self-contained cellular modem that links to our central server. Using our secured website users and caregivers can easily program their unit's medication schedule remotely, set up preferred notifications and review reports.

## Medical Alert

You can receive your MedMinder with or without medical alert. Our optional emergency alert system is embedded directly into the device and communicates with pendants worn as a necklace or watch around the home. With the push of a button the pill dispenser will open a two-way voice channel with the professionals at the UL certified monitoring center.

The wireless pendants are waterproof, feature 800 feet of transmission range, a long-life battery, and have a necklace and a male or female watch band design.

MedMinder

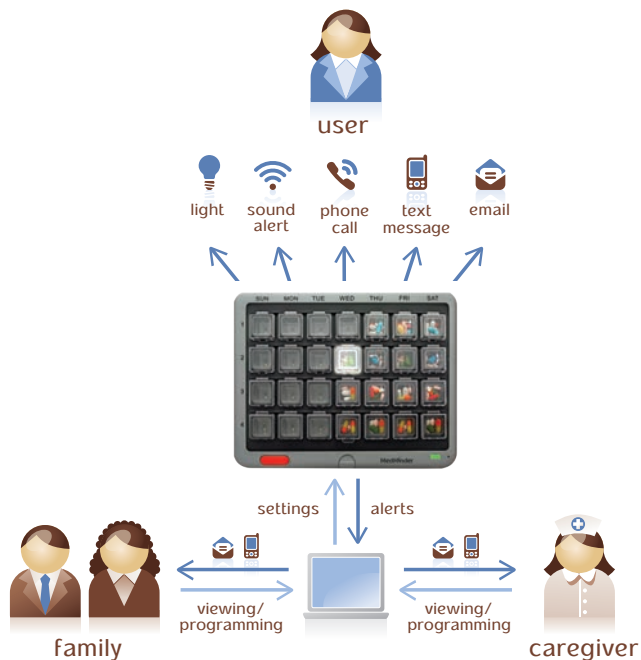

## (Re)filling

MedMinder offers two easy options to fill the medication compartments:

1. **Manually:** Simply lifting the lockable main lid gives access to all 28 compartments.
2. **Pre-fill tray:** MedMinder supplied pre-fill trays can be filled by caregiver, family or MedMinder pharmacists. Once filled they can be placed directly into MedMinder.

## Locking

MedMinder offers two types of pill dispensers: locked and un-locked. The locked pill dispenser is for those at greater risk for taking the wrong medication at the wrong time. All compartments are locked and only the specified compartment will unlock at the right time. The Jon features an “un-locking / re-locking schedule” where the dosage compartments can be setup to un-lock or re-lock before or after the scheduled dosage to allow more flexibility for the user.

## Custom Greetings

Caregivers can record voice messages at the MedMinder portal and our devices will play them so the user will hear a familiar voice. Custom greetings can be set as medication reminders, positive feedback, or even one-time personalized messages for special occasions. Caregivers can also type a personalized message that MedMinder converts to voice format.

“Don’t forget your 2:30 doctor appointment!”

“Grandpa, don’t forget to take your pills.”

“Great job, Mom! Hope you have a good day”

MedMinder™

Phone: 1-888-MED-MIND (1-888-633-6463) [www.MedMinder.com](http://www.MedMinder.com)